

Goodpac™ Technology

DISPENSE ONE GLOVE AT A TIME, AS INTENDED


When gloves are not stacked properly, dispensing becomes a problem. Very often, uneven glove packing makes gloves pulling harder, and also results in multiple glove dispensing which could be hidden cost to you especially when it ends up on the floor. Hence, a nice packaging on the outside, may not be what it seems on the inside.

Goodpac™ Technology helps to prevent uneven glove layering, and a messy packed glove box. It is a smart packing technology where every glove is alternately layered perfectly, allowing you to dispense one glove at a time, as intended.*

*All care are taken to ensure alternating layering for every single glove. However, due to our process which requires gloves to be inspected offline, some batches may require manual repacking.

Why multiple dispensing occurs in conventional packing method



Due to the messy packing method, fingers or gloves from 2nd or 3rd layer (or any other bottom layers), can overlap and be exposed at the top layer for the wearer to pull. When wearer unknowingly pulls the bottom layers, all other top layers will be pulled and dispensed as well. Hence, multiple dispensing occurs.

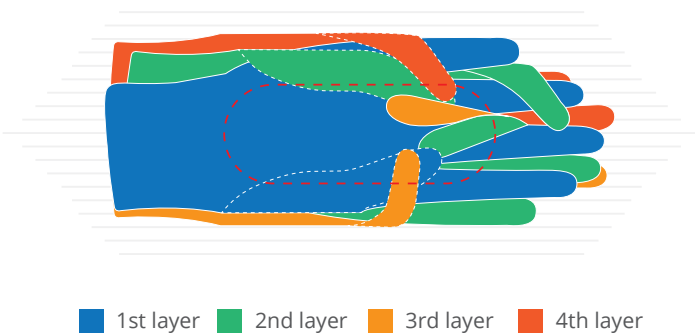

When packing, overlap of the gloves' thumbs and fingers occur as gloves are folded to fit the box


How Goodpac™ Technology dispense one glove at a time

With **Goodpac™** alternating layering method, each top layer glove cuff will always completely cover the fingers of its immediate bottom layer glove. Hence, exposing only the top most layer for the wearer to pull and therefore prevents accidentally pulling any of the bottom layered gloves.

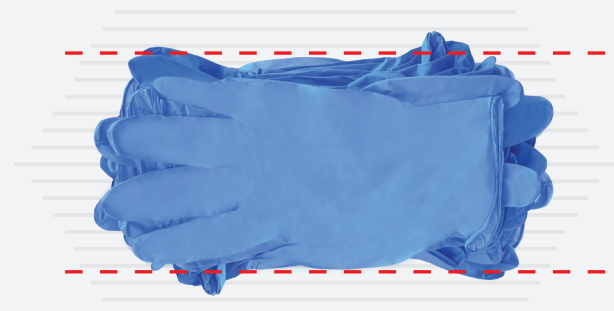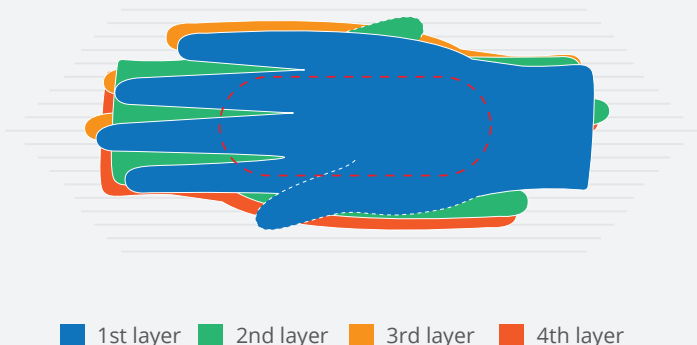

When packing, no overlap occurs as gloves are neatly stacked to fit the box

Why uneven, bulging box occurs in conventional packing method

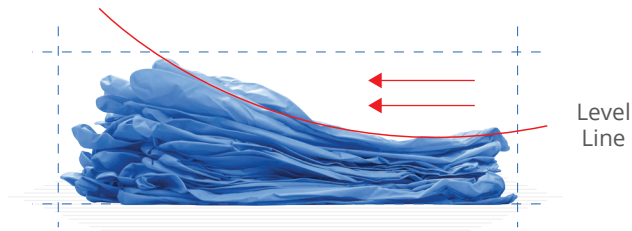


Gloves are stacked in alternate directions by clusters resulting in uneven glove layering



Bulging and clumsy looking box

How Goodpac™ Technology ensures an evenly, level surface box everytime

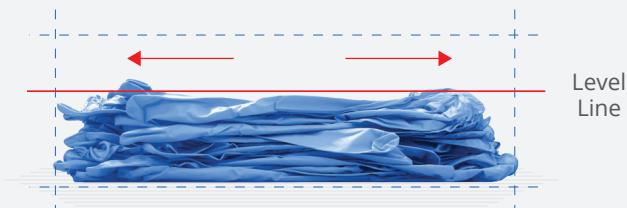

Every single glove is layered in alternate directions to provide an even glove layering



Level surface box allows stability when boxes are stacked

Benefits of Goodpac™ Technology

Enables nicely and neatly stacked gloves

Saves cost as it prevents glove wastage from multiple dispensing

Makes glove dispensing easier

Eliminates a bulging and clumsy looking glove box

Choose Goodpac™ Technology Today!



CONTACT US:

<https://hartalega.com.my/contact-us/sales-contact/>